




A Closer Look

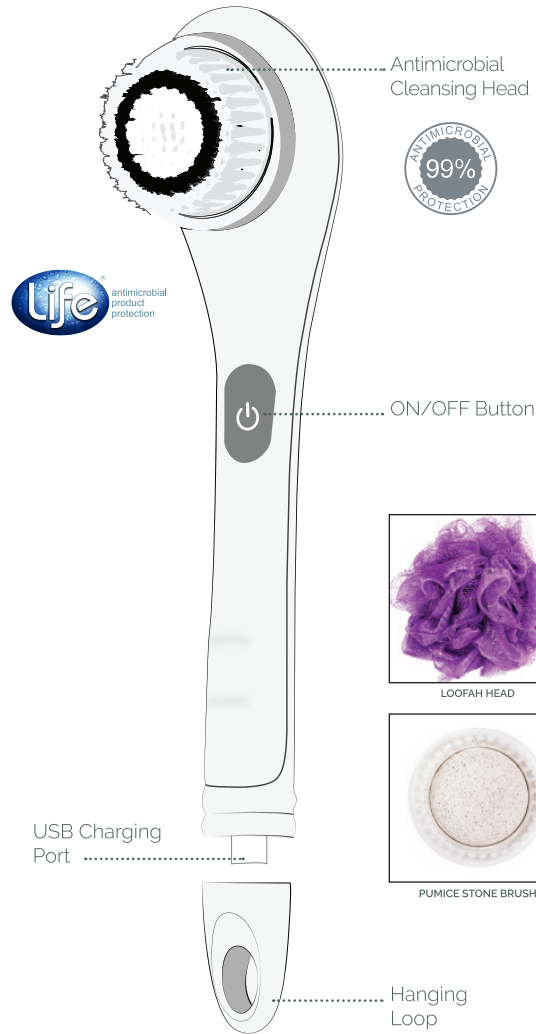

Charging

BEFORE FIRST USE




NERA is cordless and rechargeable. Before use, charge fully for 12 hours, even if the indicator light has changed color.

CHARGING

To access the USB charging port, rotate the battery cover located at the end of the NERA handle. Pull the cover out (diagram 1). Then connect the USB end of the charging cord to a USB port of a computer or a USB wall adapter and then plug into an electrical outlet. As the NERA charges, a tiny red light will come on to indicate charging is in progress, and it will turn green once your device is fully charged. Please be sure to allow the NERA to charge for the full 12 hours, even if the charging light comes off before the 12 hours, allow it to charge for the full 12 hours. Recharge the NERA only as needed and up to 12 hours only. DO NOT overcharge which can shorten the battery life. Disconnect the Charging Cord when not charging the device. NERA should be powered OFF and cannot be used while charging. Do not charge NERA for more than 24 hours.

POWERING OFF AND ON

Press the ON/OFF button to turn the device on for low speed.
 Press a second time for faster speed.
 Press a third time to turn it off.

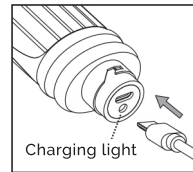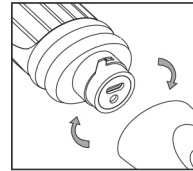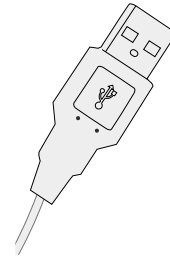


Usage Guide


READY TO USE:

Once you have assured of the full charge of your NERA device you can proceed with your treatment.


1. NERA can be used in the shower.

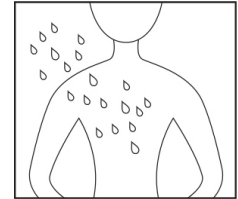

2. Apply your favorite body wash or soap.

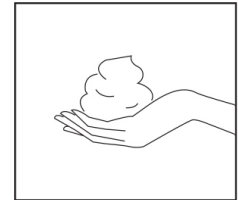


3. Place NERA where you would like to cleanse and turn on the power button and select your preferred speed out of the 2 speeds.

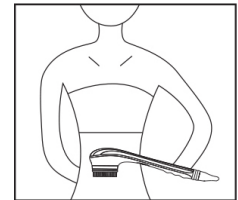


Usage Guide


4. Repeat Step 3 in all desired areas.

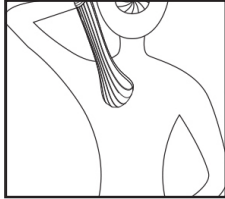

5. Rise off.

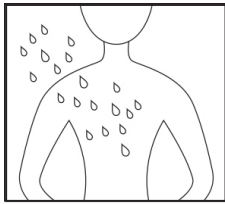

6. Clean the bath brush. Put it in a cool & dry place.

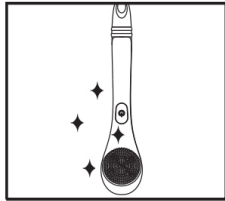

Expert Advice


SOME TIPS TO MAINTAIN YOUR NERA:

- Do not clean NERA with harsh chemical agents, and store it away from children's reach.
- Do not press down too hard as doing so may cause skin trouble or discomfort.
- Do not expose or leave NERA in direct sunlight, high heat, or boiling water, doing so may cause damage.
- For hygiene reasons, we do not recommend sharing NERA with other people.
- Exercise caution when using NERA around the eyes, ensuring that it does not touch the eyes.
- Remember to turn off the device after use. Store NERA in a clean, dry place.

BRUSH HEAD REPLACEMENT:

Put the replacement brush head on according to arrow direction.

Press in place to finish the assembling.

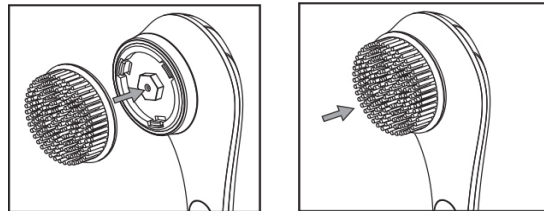

Warnings & Warranty

BASIC SAFETY PRECAUTIONS.

When using electrical products, especially when children are present, basic safety precautions should always be taken. IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF, YOU ARE WRONG. If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

DANGERS.

To reduce the risk of electrocution: DO NOT get the Product or Charging Cable wet or immerse in water. Keep this product away from open flame or flammable aerosol products. The Device is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

WARNINGS.

To reduce the risk of burns, electrocution, fire or injury to persons: Close supervision is necessary when this Product is used by, on, or near children or physically challenged persons. Use this product only for its intended use as described in this User Guide. DO NOT use this product where oxygen is being administered. The product has no serviceable parts. Never attempt to open the Device Handle for any reason. DO NOT disassemble, heat above 600°C (1400°F) or incinerate.

IMPORTANT CAUTIONS.

This is not a medical Device and is not intended to treat medical conditions. If you are under a doctor's care for any skin condition, if you think you may have a skin condition or if you have very sensitive skin, consult your physician before using this product. This product uses magnets which may affect internal medical devices such as pacemakers. Keep away from internal medical devices. This Device is electrical and contains small parts. Keep out of reach of children. THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER.

If you have any medical concerns about using NERA, please consult your dermatologist before use.

If any of the components are damaged (NERA Device Handle, or USB Charging Cord) DO NOT use. DO NOT attempt to open or services components. Email Spa Sciences customer services at customerservice@spasciences.com

Email customer service at
customerservice@spasciences.com

© 2021 Spa Sciences. All rights reserved.

Warnings & Warranty

Warranty:

Your warranty is not effective unless you register the Product within 90 days of purchase. To register go online to www.warrantyspasciences.com and complete the information required for registration.

For a period of TWO (2) years from the date of purchase, Spa Sciences LP (the "Company") warrants the NERA (the "Product") against damage, defects or failure resulting from ordinary consumer use under normal home conditions according to the instructions included with the Product (a "Covered Defect"). This warranty DOES NOT include damage, defects or failure caused by misuse or improper care of the Product; impact or insertion of foreign objects into or on the Product; application of harmful chemicals or abrasives to the Product; failure to follow manufacturer's instructions; normal wear and tear; alteration or overuse of the Product; or any professional or commercial use of the Product. If there is a Covered Defect during the two (2) years warranty period the Company agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of this limited warranty. Any claim for breach of this limited warranty must be commenced within two (2) years following the original purchase. This limited warranty applies only to a purchase of the Product from or an authorized reseller, is valid for the original end user of the Product only and cannot be transferred under any circumstances.

THE PROVISIONS OF THE FOREGOING WARRANTIES ARE IN LIEU OF ANY OTHER WARRANTY, WHETHER EXPRESS OR IMPLIED, WRITTEN OR ORAL (INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE). IN NO EVENT SHALL THE COMPANY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON ANY LEGAL THEORY. SOME STATES DO NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF SUCH IMPLIED WARRANTIES, ARE EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH HEREIN.

Additional Info

The entire operating manual including safety, care and warranty information is available at
www.neramanual.com